

ieMR Advanced

## Anaesthesia assistant log in and out

### SAAnaesthesia

#### Quick reference guide

#### How to log in to an SAA Monitor


1. While performing the Level 2 Anaesthesia check ensure data cables are connected correctly (Phillips – *Purple*, SAA – *Blue*).



2. The image below should display when you look at the SAA monitor (which will be powered up whenever the Anaesthetic machine it is connected to is powered up).



3. You should confirm your Phillips Monitor is showing your correct location e.g. induction bay 1 should display (if not, call a super user).

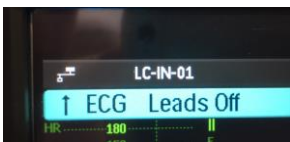

4. Use the scanner attached to the SAA monitor to scan the barcode on the wall or pendant (once) to login to the above screen / confirm the location of the SAA monitor and Phillips monitor.




wall

or



pendant

5. This should display a screen showing the correct location (which should match the Phillips monitor location) you are in/have scanned.

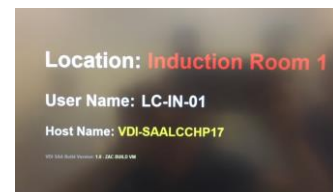

The above should lead you to a desktop with a *Purple* ieMR icon.



If this does not display, call a super user for help.

## Logging off at the end of the list/day

**Important:** At the end of the list / day / if you are moving that SAA monitor to a new location (whether it is attached to an anaesthetic machine or a free standing "bull") please close all windows (SAA and Powerchart) and double click on the "Logout" icon on the desktop. You do not need to power down the SAA monitor, it remains on.

